

INSTRUCTION FOR USE

SPMM/RA/TF /05/INS/IFU, Rev. 06 Date 11/05/2026



INSTRUCTION FOR USE

Insulin Syringe

		
DATE: 11/05/2026	DATE: 11/05/2026	DATE: 11/05/2026
Prepared By Executive- QA	Reviewed By RA Head	Approved By Manging Director/ Authorized signatory

INSTRUCTION FOR USE

SPMM/RA/TF /05/INS/IFU, Rev. 06 Date 11/05/2026

Follow manufacturer's instructions for use and institutional procedures for administration of Insulin

PRODUCT NAME: Insulin Syringe

DESCRIPTION

The SPM's Insulin Syringe is a plastic syringe designed for subcutaneous injection of a desired dose of insulin. The insulin syringe consists of a graduated barrel, plunger rod and needle/hub assembly. The insulin syringes are sterile, single use, and non-pyrogenic. This device operates on the same principle as piston syringes. It is sterilized by ethylene oxide (EO).

Basic UDI-DI: 89043798ISHP

PRODUCT IMAGE



DEVICE VARIANTS (Models)

Product Description	Gauge Size (Needle Size)
Insulin Syringe U-40	29G × 6 mm
Insulin Syringe U-40	30G × 8 mm
Insulin Syringe U-40	31G × 12.7 mm
Insulin Syringe U-100	29G × 6 mm
Insulin Syringe U-100	30G × 8 mm
Insulin Syringe U-100	31G × 12.7 mm
Insulin Syringe U-40	29G × 6 mm
Insulin Syringe U-40	30G × 8 mm
Insulin Syringe U-40	31G × 12.7 mm
Insulin Syringe U-100	29G × 6 mm
Insulin Syringe U-100	30G × 8 mm
Insulin Syringe U-100	31G × 12.7 mm

Available in red and orange colour.

INSTRUCTION FOR USE

SPMM/RA/TF /05/INS/IFU, Rev. 06 Date 11/05/2026

INTENDED USE (Purpose)

The Insulin Syringe is intended for subcutaneous administration of insulin in diabetic patients for transient and single use only in the insulin-dependent Diabetic patients.

INDICATIONS FOR USE

The insulin syringe is intended for use in diabetic patients for the manual aspiration and subcutaneous injection of insulin and is intended for single use only.

INTENDED PATIENT POPULATION

The Insulin Syringe is intended for diabetic patients irrespective of age and gender.

INTENDED USER

This device is intended for use by healthcare professionals. Self-administration by trained patients is permitted only when instructed by a qualified healthcare professional

INTENDED USE ENVIRONMENT

Home setting, hospital, or clinical setting.

STERILISATION

The device is supplied sterile and non-pyrogenic, sterilised using ethylene oxide (EO).

CLINICAL BENEFITS

Short and thin attached needle for more safety and comfort.
Low dead space for more accurate dosage and minimal wastage.
Bold and clear graduation markings for ease of reading.

PRINCIPLE OF OPERATION

The SPM's insulin syringe is used for the administration of insulin via a fine needle. It consists of a marked barrel, plunger, and fine needle with a protective needle cap. The user draws insulin into the syringe, removes air bubbles, and injects it into the fatty tissue under the skin to enable absorption into the bloodstream and regulation of blood glucose levels.

INSTRUCTION FOR USE

SPMM/RA/TF /05/INS/IFU, Rev. 06 Date 11/05/2026

PERFORMANCE CHARACTERISTICS

Provides safe and accurate insulin administration for diabetic patients.
Low dead space design for accurate dosing and minimal drug wastage.
Available in red and orange colour options.
Bottom cap and needle cap enhance handling and safety.

DEVICE COMPONENT DESCRIPTION (Materials)

Barrel – Polypropylene (PP)
Plunger – Polypropylene (PP)
Gasket – Thermoplastic Elastomer (TPE)
Needle Cap – Polyethylene (PE)
Plunger (Bottom) Cap – Polyethylene (PE)
Lubricant – Silicone oil
Adhesive – UV glue
Needle/Hub Assembly – Stainless steel cannula with integral hub

COMPONENTS OF THE DEVICE

Barrel (graduated)
Plunger rod
Gasket
Needle with hub (attached)
Needle cap
Plunger bottom cap

CONTRAINDICATIONS

Sr. No	Contraindications	Description
1	Allergy to materials used in syringe	Patients with known allergies or sensitivities to syringe materials (e.g., specific plastics, latex, or needle metals) should avoid use or take precautions.
2	Infection	In patients with active infections or an increased infection risk, using a contaminated or improperly sterile syringe could worsen the condition.
3	Coagulopathy or bleeding disorders	Patients with bleeding disorders may have a higher risk of bleeding or haematoma formation at the injection site.
4	Existing vascular conditions	Patients with vascular or anatomical abnormalities may face increased risks of complications or require specialized equipment.
5	Incompatible medications or fluids	These syringes are designed only for insulin ; using incompatible drugs may reduce efficacy or trigger adverse reactions.
6	Device misuse	Damaged or incorrectly handled syringes can lead to improper dosing, tissue injury, or insulin degradation; use strictly as per instructions.

INSTRUCTION FOR USE


SPMM/RA/TF /05/INS/IFU, Rev. 06 Date 11/05/2026

POTENTIAL COMPLICATIONS

- Hypoglycaemia due to overdosing or incorrect administration.
- Local pain, bruising, or bleeding at injection site.
- Infection or abscess if asepsis is not maintained.
- Lipohypertrophy/lipoatrophy with repeated injections at the same site.
- Allergic reactions (to materials or insulin).
- Needle-stick injury to user or healthcare professional.
- Cross-contamination with reuse.

DIRECTIONS FOR USE

- Verify packaging: Use only if the sterile pack is intact and within expiry.
- Prepare dose: Confirm the insulin concentration (U-40 or U-100) matches the syringe type. Wash hands.

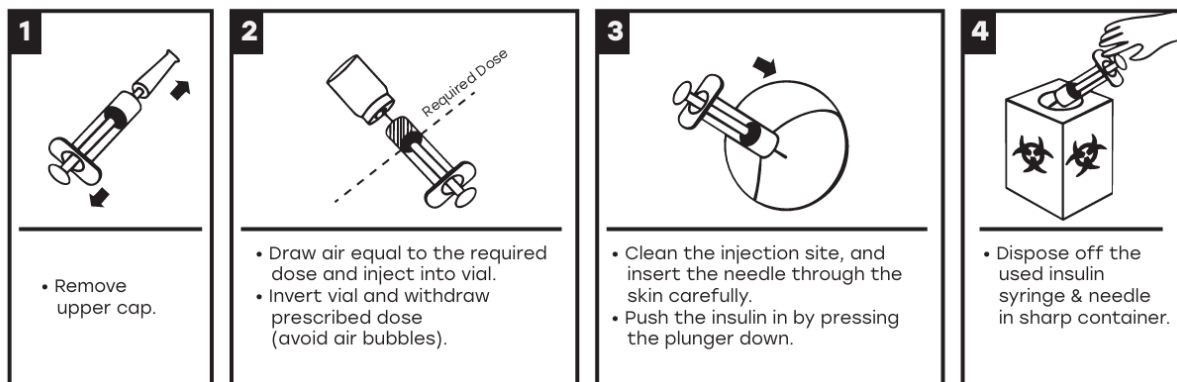

HOW SUPPLIED

The SPM's Insulin Syringe is supplied sterile, single-use, in sterile Blister plastic paper/ Ribbon / Multi-pack pouch packaging, complete with instructions for use. Refer to the device label for batch or lot identification.

WARNINGS

- Improper transport/handling may damage the device or packaging.
- Product is non-toxic, sterile & non-pyrogenic if the package is unopened and undamaged.
- Single use only. Do not re-sterilise, reuse, or re-process.
- Re-sterilising can compromise performance and increase cross-contamination risk.
- Reuse may cause cross-infection.
- Remove the needle cap carefully to avoid needle-stick injury.
- Recommended for Personal Use.

INSTRUCTION FOR USE

SPMM/RA/TF /05/INS/IFU, Rev. 06 Date 11/05/2026

- Do not inject air into the patient; avoid manipulating the plunger in a way that affects dose accuracy.
- Reuse of this device may result in loss of sterility, inaccurate dosing, infection, cross-contamination, or serious injury.

PRECAUTIONS

- Read instructions before use.
- Store in a clean, dry place; avoid direct sunlight; keep at room temperature.
- Do not use after the expiry date.
- Do not use if sterile packaging is damaged/opened prior to use.
- Use immediately after opening.
- Always match syringe type to insulin concentration (U-40 with U-40; U-100 with U-100).

CAUTIONS

- For single use only. Do not reuse.
- Do not use if packaging is damaged or opened.
- Do not use after expiry date.
- Do not re-sterilize.
- Discard immediately after use in an approved sharps container.
- Reuse may cause infection or injury.
- Consult a physician before changing syringe type or insulin mixing method.
- Keep out of reach of children.
- Store in a cool, dry place.

ADVERSE EVENTS

- Phlebitis (not typical for s.c. injection; generally, not expected)
- Infection / toxic reaction
- Injury / needle-stick injury
- Incomplete dose administration
- Blood-glucose instability (e.g., hypoglycaemia)
- Bruising/haematoma at injection site
- Lipodystrophy with poor site rotation
- Allergic reactions
- Cannulation failure (not applicable; syringe only)

STORAGE CONDITION

Store at a controlled temperature (10–40°C) in a cool, dry, dark place.

Note: Avoid exposure to direct sunlight. Exceeding the recommended storage conditions may affect product quality.

INSTRUCTION FOR USE

SPMM/RA/TF /05/INS/IFU, Rev. 06 Date 11/05/2026

SHELF LIFE

5 years from the date of manufacture. Do not use after the indicated sterility expiry date.

RETURN OF DAMAGED PRODUCT

Return the product in its original packing identified by the batch number, purchase information, your reference and reason for return. Please contact your local distributor office regarding product return/exchange.

DISPOSAL

- Dispose of the used syringe immediately after use in an approved sharps container. Do not recap the needle.
- After use, this product may be a potential biohazard.
- Dispose of after single-patient use during a single procedure.
- Dispose of the used product and packaging as biomedical waste according to hospital procedures and local regulations.
- Any open/unused items from opened packs should be disposed of as biomedical waste.
- Leftover products should be disposed of as biomedical waste.

REPORTING

Any serious incident occurring in relation to the device shall be reported to the manufacturer and to the competent authority of the Member State in which the user and/or patient is established, in accordance with Regulation (EU) 2017/745.

LEGAL MANUFACTURER AND MANUFACTURING FACILITY

<u>LEGAL MANUFACTURER AND MANUFACTURING FACILITY</u>	<u>AUTHORIZED REPRESENTATIVE</u>
 Manufacturer: SPM Medicare Pvt, Ltd B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India. Manufacturing Site 1- B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India. Manufacturing Site 2- Industrial Building Plot no 26, Sector 158, Noida, Gautam buddha, Nagar, Uttar Pradesh, 201305. Customer care: No.: +91-840-7070-718 Email: info@spmmedicare.com Web Site: www.spmmedicare.com	EC REP CMC MEDICAL DEVICES & DRUGS S.L. C/ Horacio Lengo, 18 29006, Málaga, Spain Customer Care No.: +34 951 214 054 Email: info@cmcmedicaldevices.com

INSTRUCTION FOR USE

SPMM/RA/TF /05/INS/IFU, Rev. 06 Date 11/05/2026

EXPLANATION OF SYMBOLS:

The symbols used on the device label and packaging comply with EN ISO 15223-1.

Symbol	Explanation	Symbol	Explanation
	Caution		Catalogue Number
	Manufacturer		Batch code
	Date & Country of manufacturer		Use-by date
	Do not re-sterilize		Temperature Limit
	Consult instructions for use		Do not use if package is damaged and consult <i>instructions for use</i>
	Keep away from sunlight		Non-Pyrogenic
	Do not re-use		Sterilized using EO
	Keep dry		Handle with care
	This side up		Maximum stacking
	Do not use hook		Recycling

INSTRUCTION FOR USE

SPMM/RA/TF /05/INS/IFU, Rev. 06 Date 11/05/2026

	Medical Device		Model number
	Unique Device Identifier		Single Sterile Barrier System
	Single Sterile Barrier System with protective packaging outside		Biological risks

Note:

- * Not made with natural rubber latex.
- * Not made with DEHP.

Revision History

Rev. No	Section updated/Changed	Description of Change	Date	Approval
05	Cautions	Added a caution regarding single use and updated the wording in accordance with compliance requirements.	14/04/2026	Managing Director/Auth.signatory
06	LEGAL MANUFACTURER AND MANUFACTURING FACILITY	Change of Authorized Representative.	11/05/2026	Managing Director/Auth.signatory