

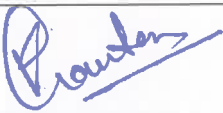
INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ MV Set/ 04. - Rev. 05 Date: 11-05-2026



INSRUCTION FOR USE

MEASURED VOLUME SET (BURETTE SET)

 DATE: 11-05-2026	 DATE: 11-05-2026	 DATE: 11-05-2026
Prepared By Executive- QARA	Reviewed By Head-RA	Approved By Managing Director/Auth. Signatory

INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ MV Set/ 04. - Rev. 05 Date: 11-05-2026

Follow manufacturer's instructions for use and institutional procedures of Closed IV Cannula

Product Name- Measured Volume Set (Burette Set)

DESCRIPTION

The Measure Volume Set (Burette Set) is used to administer Intravenous fluid and medicines into human circulating system by using intravenous catheter or cannula. The Measure Volume Set (Burette Set) is specially designed to administer measured volume of infusion fluid through gravity as well as pressure pump method in children. The product is sterilized using EO (Ethylene Oxide) gas. The use of the product is restricted to a qualified doctor or a paramedic.

IMAGE OF THE DEVICE



DEVICE VARIANTS:

Sr. No.	Model	Product Code
1.	ACCU PEDIA MEASURED VOLUME IV SET REGULAR - 110ML, RIBBON PACK	2CRB2
2.	ACCU PEDIA MEASURED VOLUME IV SET REGULAR- 150ML, RIBBON PACK	2CRB6
3.	ACCU PEDIA MEASURED VOLUME IV SET PLUS- 110ML, PEEL OPEN PACK	2CPB1
4.	ACCU PEDIA MEASURED VOLUME IV SET PLUS - 150ML, PEEL OPEN PACK	2CPB5
5.	ACCU PEDIA MEASURED VOLUME IV SET REGULAR - 110ML, PEEL OPEN POUCH	2BKPB3
6.	ACCU PEDIA MEASURED VOLUME IV SET REGULAR - 150 ML, PEEL OPEN POUCH	2BKPB4
7.	EASY PEDIA MEASURED VOLUME IV SET - 110ML, HM PACK	2ACHB2
8.	EASY PEDIA MEASURED VOLUME IV SET - 150ML, HM PACK	2ACHB6



INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ MV Set/ 04. - Rev. 05 Date: 11-05-2026

INTENDED PURPOSE

MEASURED VOLUME SET Equipped with a micro drip with reduced drop size of 60 drop per ml and a burette type chamber of 100ml or 150ml capacity for infusion fluids intravenous to pediatric patient with the help of hypodermic device.

INDICATIONS FOR USE

- Infusion of IV Fluids
- To prevent/ treat dehydration
- Infusion of IV Drug Administration
- **INTENDED PATIENT POPULATION**
- The Measure Volume Set (Burette Set) is intended for administer intravenous fluid and/or medication to a patient intravenously.

- **INTENDED USER**

- The intended users for The Measure Volume Set (Burette Set) are healthcare professionals trained in performing Infusion.

- **INTENDED USE ENVIRONMENT**

- The Measure Volume Set (Burette Set) is intended to be used in hospitals and clinics.

- **CLINICAL BENEFITS**

- The Measure Volume Set (Burette Set) improves clinical care by providing secure and consistent intravenous access, reducing the need for frequent line replacements and minimizing infection risks. It ensures efficient fluid management and patient safety.

Intended Patient Population:

The MV Set is intended for use in neonatal, pediatric, and adult patients requiring accurate volume-controlled infusion of IV fluids or medications.

DURATION OF USE:

Short term Use (Less than 30 Days).

STORAGE CONDITIONS:

Store at specified temperature and humidity (Temperature in between 10°C to 40°C and Humidity in between 55RH to 75RH).

CONTRAINDICATIONS/WARNINGS:

The instructions should be followed carefully before use.

- DO NOT re-sterilize and /or reuse.
- Re use of single use device creates a potential risk for patient on use. It may



INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ MV Set/ 04. - Rev. 05 Date: 11-05-2026

- lead to contamination and /or impairment of functional capability.
- Do not use on injured or infected surface.
- Check integrity at the seal where top cap and barrel meet to ensure sterility
- Visually inspect and carefully check the product and packaging before use.
- Do not use with light sensitive or photo sensitive drugs & paclitaxel type of
- chemotherapeutic drugs. Keep away from children

INDICATION/INTENDED APPLICATION:

To prevent or treat dehydration.

PRINCIPLE OF OPERATION:

This is a disposable measured volume Infusion Set of different variants. This device is supplied sterile in paper pouch. Based on stability studies, the shelf life for sterile products is recommended to 5 years from the date of sterilization/manufacturing. The product life is mentioned on the product label.

PERFORMANCE CHARACTERISTICS:

The Measured Volume Set delivers a controlled micro-drip rate of 60 drops/mL for precise infusion.

It features a calibrated burette chamber (110 mL/150 mL) for accurate volume measurement. Equipped with an air vent and flow regulator to ensure safe, bubble-free infusion. The device is EO sterilized and tested for leak-proof performance and biocompatibility.

HOW SUPPLIED:

MV Set is supplied sterile, Single unit is packed in individual Ribbon/HM Pouch made of Low Density Polyethylene (LDPE) as primary packaging.

-then 10 unit of single packed pouches are in a master low density polyethylene pouch as a secondary packaging.

-Then 08 units of secondary packed pouches (80 units of the product) are supplies in a 5 Ply Corrugated Box as a tertiary packaging with external identification stickers.

-Product IFU listed on the primary pack.

MATERIALS USED:

S. NO.	Components	Material
1	Spike Cap	Poly Propylene
2	Spike	ABS
3	Drip Chamber	PVC-DEHP Free
4	Button Filter	PP
5	P.V.C. Tube	P.V.C.-DEHP Free
6	Latex Bulb	Latex
7	Roller Body	Polypropylene
8	Roller	ABS
9	Burette/ Cylinder	Poly Propylene



INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ MV Set/ 04. - Rev. 05 Date: 11-05-2026

10	Vent Cap	Poly Propylene
11	Pinch Cock	Polypropylene
12	Stopper	Silicon
13	Micro Pin	SS

POTENTIAL COMPLICATIONS:

- Air embolism due to improper priming or air bubble entry.
- Infection risk if aseptic technique is not followed.
- Inaccurate fluid administration from malfunction or flow regulator error.
- Allergic reaction to material components (e.g., latex).



CAUTIONS:

- Store in a cool and dry place.
- Discard after single use.
- Do not use if package is open or damaged.
- Do not store at extreme temperature and humidity.
- Do not use with blood and blood components.
- Do not use device after shelf life expires.

PRECAUTIONS

- Do not use if packing is damaged.
- Do not use if protective caps are loose or missing.
- Use injection site only for injecting medicine.
- Avoid introduction of air bubbles to prevent air embolism.

INSTRUCTIONS FOR USE

- Check the packing carefully, if packing is found damaged, torn or pierced, discard the piece.
- Wash-up and scrub hands and preferably use pre-sterile protective gloves.
- Peel opens the pack and draw out the device aseptically.
- Close the flow regulator.
- Remove protective cap from plastic spike & insert the spike at the top of the bottle to its full length.
- Squeeze the drip chamber until it is half filled.
- Removing prime stop connector from rotating luer locks connector.
- Open flow regulator & allow solution to pass until all air bubbles in the tube are removed. Close the flow regulator.
- Connect the set to intended device and regulate flow by using flow regulator.
- Dispose off the used device as per local applicable laws.

SHELF LIFE

The expiration date of the product is the sterility expiration date. The product should not be use after the indicated sterility expiration date. Recommended shelf life of the product is 05 years from the date of manufacturing.

INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ MV Set/ 04. - Rev. 05 Date: 11-05-2026

RETURN OF DAMAGED PRODUCT:

Return the product in its original packing identified by the batch number, purchase information, your reference and reason for return. Please contact your local distributor office regarding product return/exchange.

DISPOSAL OF THE USED DEVICE

- Disposal of the used product should be done according to hospital procedures, State Regulations, national legislation or regional legislation.
- Device should be disposed after use on individual patient during a single procedure.
- Disposal of the packaging material should be done as biomedical waste.
- Any open packaging should be disposed as biomedical waste.
- Leftover products should be disposed as biomedical waste

REPORTING

In case of any adverse events and/or potentially sight-threatening complications that may reasonably be regarded as product-related and that were not previously expected in nature, severity or occurrence, contact the manufacturer or authorized representative and competent authority of the member state where the user and/or patient is established.

<u>LEGAL MANUFACTURER AND MANUFACTURING FACILITY</u>	<u>AUTHORIZED REPRESENTATIVE</u>
 Manufacturer: SPM Medicare Pvt, Ltd B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India. Manufacturing Site 1- B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India. Manufacturing Site 2- Industrial Building Plot no 26, Sector 158, Noida, Gautam buddha, Nagar, Uttar Pradesh, 201305. Customer care: No.: +91-840-7070-718 Email: info@spmmedicare.com Web Site: www.spmmedicare.com	EC REF CMC MEDICAL DEVICES & DRUGS S.L. C/ Horacio Lengo, 18 29006, Málaga, Spain Customer Care No.: +34 951 214 054 Email: info@cmcmmedicaldevices.com

EXPLICATION OF SYMBOLS

Symbol	Meaning	Symb ol	Meaning
	Caution/Warning (Read instruction before use)		Catalogue Number

INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ MV Set/ 04. - Rev. 05 **Date:** 11-05-2026

	Manufacturer Logo		Batch/Lot Number
	Date and country of manufacturer		Date of Expiry
	Do not re-sterilize		Temperature Limitation
	Humidity Limitation		Do not use if package damaged
	Consult IFU		Non Pyrogenic
	Keep away from sunlight		Sterilized using EO
	Don't re-use		Handle with care
	This side up		Maximum stacking
	Do not use hook		Conformity of European Norm with Notified Body Number
	Medical Device		Model Number
	Unique Device Identifier		Single Sterile Barrier System
	Single Sterile Barrier System with protective packaging outside		Authorized European Representative
	Keep Dry		



INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ MV Set/ 04. - Rev. 05 **Date:** 11-05-2026

Revision History:

Rev. No	Reason for Change	Date	Approval
05	Change of Authorized Representative.	11/05/2026	Managing Director