

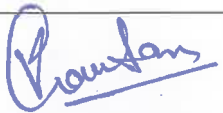
INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ 3 Way Stop Cock/ 06. - Rev. 04 **Date:** 11-05-2026



INSTRUCTION FOR USE

Extension line with Three-way stopcock

 DATE: 11-05-2026	 DATE: 11-05-2026	 DATE: 11-05-2026
Prepared By Executive- QARA	Reviewed By Head-RA	Approved By Managing Director/Auth. Signatory

Follow manufacturer's instructions for use and institutional procedures for Extension line with Three-way stopcock

PRODUCT NAME: Extension line with Three-way stopcock

DESCRIPTION:

The SPM's Extension line with Three-way Stopcock is a sterile, single-use, Pyrogen-free Extension line with integrated 3-way stopcock. The Extension line is available in different lengths ranging from 10cm to 250cm PVC Tube. The Extension line consists of an integrated 3-Stopcock at one end and a male luer at another end. The Stopcock contains Rotating knob and two female ports which allow multiple drug administration. The device is ethylene oxide sterilized.

IMAGE OF THE DEVICE:



DEVICE VARIANTS:

Sr. No.	Model	Product Code
1.	ACCU LINE EXTENSION LINE - 50CM, PEEL OPEN POUCH	2FPPA
2.	ACCU LINE EXTENSION LINE - 100CM, PEEL OPEN POUCH	2FPPB
3.	ACCU LINE EXTENSION LINE - 150CM, PEEL OPEN POUCH	2FPPC
4.	ACCU LINE EXTENSION LINE - 200CM PEEL OPEN POUCH	2FPPD
5.	NEEDLE FREE VALVE WITH EXTENSION LINE ONE WAY	2GPNU
6.	NEEDLE FREE VALVE WITH EXTENSION LINE ONE WAY	2GPND
7.	NEEDLE FREE VALVE WITH EXTENSION LINE ONE WAY	2GPNT

INTENDED PURPOSE:

The Extension line with 3-way stop cock is intended for infusion therapy. Extension line is used to extend the infusion line during administration of IV fluid. It also provides an alternative channel for the introduction of medication to the patient.

INDICATIONS FOR USE:

- SPM's Extension line with three-way stopcock allows connection of multiple IV lines to a single access point, enabling fluid and medication administration without additional venepunctures.
- It is often used for invasive pressure monitoring in critical care settings, allowing for the monitoring of blood pressure or central venous pressure through one of the ports.

INTENDED PATIENT POPULATIONS:

The extension line 3-ways stop cock is intended for children, Neonatal, adolescents, and adults who are healthy and seeking IV Therapy Patients, Dialysis, needing regular infusions, connection of multiple IV lines or syringes to a single access point, enabling the administration of different medications or fluids without needing to insert multiple needles.

INTENDED USERS:

Extension line with three-way stopcock is intended to be used by the trained healthcare professionals including doctors, registered practitioners, and paramedics or nursing staff.

INTENDED USE ENVIRONMENT:

Extension line with three-way stopcock is intended to use in hospitals and clinics

CLINICAL BENEFITS:

- The three arm handle with arrow indication for smooth and easy operation
- 6% taper in luer lock allows connection with a variety of medical devices.
- Clear Luer lock ABS material and kink-resistant PVC (Latex-free) extension line.
- Allows Multiple drugs/ IV fluid administration

PRINCIPLE OF OPERATION:

The Operation of Extension line with three-way stopcock is relies on its three-port valve system, which allows for the control of fluid flow. The stopcock can be rotated to open or close specific ports, enabling multiple dose administration.

PERFORMANCE CHARACTERSTICS:

- Interoperability with other vascular access devices.
- Allow multiple dose administration

MATERIALS USED:

- Extension Tube - PVC
- Stopcock Body - Poly Carbonate
- Female Ports - Poly Propylene
- Rotating Knob - Poly Propylene
- Luer Connector- HDPE
- Protecting cap- Poly Propylene

CONTRAINDICATIONS:

Sr. No	Contraindications	Description
1	Allergy to Extension line with 3-way Stopcock Materials	Patients with known allergies or sensitivities to materials in the Extension line with 3-ways stopcock (e.g., specific plastics) should avoid using the set to prevent allergic reactions.
2	Infection	In patients with active infections or conditions that increase infection risk, using a contaminated or improperly sterile Extension line with 3-way stopcock could exacerbate their condition.
3	Incompatible Medications or Fluids:	Avoid using the Extension line with 3-ways stop cock for fluids or medications that are incompatible with the materials of the set, as this could lead to chemical reactions or reduced efficacy.
4	Existing Vascular Conditions:	Patients with specific vascular conditions or anatomical abnormalities may be at increased risk of complications or may require different Extension line with 3-ways stop cock equipment.
5	Device Misuse:	If the Extension line with 3-ways stop cock is damaged or used incorrectly, it could lead to improper dosing or injury. Ensure the set is used according to manufacturer instructions.

POTENTIAL COMPLICATIONS:

Fluid leakage and cross-contamination, Difficulty turning, Incompatibility with medications, Blockage and obstruction, Sterility breakdown, Misalignment of ports, Material deterioration, Limited durability in high-pressure.



CAUTIONS:

- Store in a cool and dry place
- Discard after single uses
- Do not use if the package is open or damage
- Do not store at extreme temperature and humidity
- Dispose the blood contaminated product as per local applicable laws.

EXPIRATION PERIOD: The device has a shelf life of 5 years from the date of manufacturing. The product should not be used after the indicated sterility expiration date.

STORAGE: Store at a controlled temperature (10-40°C). Do not store at freezing temperature.

INSTRUCTIONS FOR USE:

- **Step 1.** Peel Open the pouch.
- **Step 2.** Unscrew the luer cap of three way stop cock female luer.
- **Step 3.** Connect male portion of I.V. Set with female luer of Three Way Stop cock.
- **Step 4.** Remove luer cap of male luer tubing.
- **Step 5.** Hold male luer of tubing upright and open clamp of I.V. Set and allow the solution to displace air entirely from tubing.

- **Step 6.** Close clamp.
- **Step 7.** Connect male luer of tubing to venepuncture device.

HOW SUPPLIED

The Extension line with 3-way stopcock is supplied in sterile, single-use Paper pouch packaging, complete with instructions for use and necessary accessories (if any).

WARNINGS

- **For Professional Use Only:** This device must be used only by trained healthcare professionals
- **Sterile Packaging:** Do not use if the packaging is damaged or opened prior to use.
- **Single Use Only:** The Extension line with 3 way stop cock is designed for single use. Do not reuse or re-sterilize.
- **Proper Handling:** Handle with care to avoid damage to the device, which could compromise its sterility and functionality
- **Allergic Reactions:** Check patient history for allergies to materials in the Extension line with 3way stopcock.
- **Inspect Before Use:** Ensure all components are present and undamaged before use
- **Follow Instructions:** Carefully follow the included instructions for use to avoid complications
- **Dispose Properly:** Dispose of the device in accordance with local regulations for medical waste

PRECAUTIONS:

- Inspect the stopcock and extension line for cracks, leaks, or signs of damage before use.
- Ensure that all luer connections are secure and leak-proof before beginning infusion.
- Avoid applying excessive torque while turning the stopcock handle.
- Do not force the stopcock handle beyond its intended rotation limits to prevent breakage.
- Avoid over-pressurizing the device beyond its specified pressure rating.
- Replace immediately if any malfunction or obstruction is observed during use.
- Use aseptic technique while connecting/disconnecting to minimize infection risk.
- Flush the extension line before use to remove air and avoid air embolism.

RETURN OF DAMAGED PRODUCT:

Return the product in its original packing identified by the batch number, purchase information, your reference and reason for return. Please contact your local distributor office regarding product return/exchange.

DISPOSAL OF THE USED DEVICE:

- Disposal of the used product should be done according to hospital procedures, State Regulations, National legislation or regional legislation.

- Device should be disposed after use on individual patient during a single procedure.
- Disposal of the packaging material should be done as biomedical waste.
- Any open packaging should be done as biomedical waste.
- Leftover products should be disposed as biomedical waste.

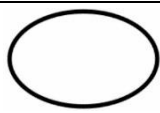
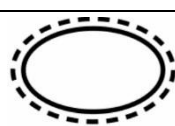
REPORTING:

In case of any adverse events and/or potential threat/ complications that may reasonably be regarded as product related and that were not previously expected in nature, severity or occurrence, contact the manufacturer or authorized representative and competent authority of the member state where the user and/or patient is established.

LEGAL MANUFACTURER AND MANUFACTURING FACILITY	AUTHORIZED REPRESENTATIVE
 Manufacturer: SPM Medicare Pvt, Ltd B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India. Manufacturing Site 1- B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India. Manufacturing Site 2- Industrial Building Plot no 26, Sector 158, Noida, Gautam buddha, Nagar, Uttar Pradesh, 201305. Customer care: No.: +91-840-7070-718 Email: info@spmmedicare.com Web Site: www.spmmedicare.com	EC REP CMC MEDICAL DEVICES & DRUGS S.L. C/ Horacio Lengo, 18 29006, Málaga, Spain Customer Care No.: +34 951 214 054 Email: info@cmcmmedicaldevices.com

EXPLANATION OF SYMBOLS:

Symbol	Explanation	Symbol	Explanation
	Caution		Catalogue Number
	Manufacturer		Batch code
	Date & Country of manufacturer		Use-by date
	Do not re-sterilize		Temperature Limit

	Consult instructions for use		Do not use if package is damaged and consult <i>instructions for use</i>
	Keep away from sunlight		Non-Pyrogenic
	Do not re-use		Sterilized using EO
	Keep dry		Handle with care
	This side up		Maximum stacking
	Do not use hook		Authorized representative in the European Community
	Medical Device		Model number
	Unique Device Identifier		Single Sterile Barrier System
	Single Sterile Barrier System with protective packaging outside		Biological risks

Revision History:

Rev. No	Reason for Change	Date	Approval
04	Change of Authorized Representative.	11/05/2026	Managing Director