



INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ PFS with sterile water/ 46 - Rev. 02 Date:11-05-2026



INSTRUCTION FOR USE

Pre filled syringe with Sterile water

 DATE: 11-05-2026	 DATE: 11-05-2026	 DATE: 11-05-2026
Prepared By Executive- QARA	Reviewed By Head-RA	Approved By Managing Director/Auth. Signatory

INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ PFS with sterile water/ 46 - Rev. 02 Date:11-05-2026

PRODUCT NAME: Prefilled syringe with Sterile water

Follow manufacturer's instructions for use and institutional procedures for Foley Catheter Balloon Inflation using Pre-filled syringe with sterile water.

Description:

Pre-filled Syringe with Sterile water is a sterile, clear, polypropylene luer lock syringe pre-filled with sterile water. It is intended to inflate the indwelling Foley catheter balloon after insertion, securing it in place. The syringe contains sterile water, ensuring contamination-free use and is intended for single-use only. The product is sterilized by Gamma irradiation.

IMAGE OF THE DEVICE



Fig. Pre-filled syringe with sterile water

MODELS:

Sr. No.	Device / Model / Size / Variant	Catalogue No. / UDI-DI / Remarks
1.	Pre- filled syringe with sterile water, 5ml	1DR7B5S
2.	Pre- filled syringe with sterile water, 10ml	1DR7BXS
3.	Pre- filled syringe with sterile water, 20ml	1DR7PYS
4.	Pre- filled syringe with sterile water, 30ml	1DR7PCS
5.	Pre- filled syringe with sterile water, 50ml	1DR7PZS

INTENDED USE:

The Pre-filled Syringe with sterile water is intended for the safe and efficient inflation of indwelling urinary Foley catheter balloon, to retain the catheter in the bladder.

INDICATIONS FOR USE

The Pre-filled Syringe with Sterile Water is indicated for inflating the retention balloon of Foley catheters during urinary catheterization procedures in a sterile clinical setting.

INTENDED USE ENVIRONMENT

The device is intended to be used in hospitals, clinics, emergency departments, long-term care facilities, and outpatient surgical centers.

CLINICAL BENEFITS

- Ensures accurate and safe balloon inflation

INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ PFS with sterile water/ 46 - Rev. 02 Date:11-05-2026

- Eliminates the need for manual filling
- Maintains sterility and reduces risk of contamination
- Saves time and improves procedural efficiency

PRINCIPLE OF OPERATION

The device operates through manual actuation of the plunger, which dispenses pre-measured sterile water via a luer-lock connection to inflate the catheter balloon, anchoring it in the bladder.

PERFORMANCE CHARACTERISTICS

- Provides accurate, pre-measured volumes (e.g., 5 mL, 10 mL, etc.) for balloon inflation.
- Luer lock design ensures secure connection with Foley catheter ports.
- Sterile water remains contamination-free due to gamma sterilization.
- Smooth plunger movement supports controlled and consistent delivery of contents.

Material Used:

Raw Material	Component / Use
Polypropylene	Barrel
TPE (Thermoplastic Elastomer)	Gasket
Silicon 1000 CST	Gasket Mixing
Polypropylene	Plunger rod
Polypropylene	Luer Cap
Sterile Water	Saline Solution

CONTRAINDICATIONS:

- The Pre-filled Syringe with Sterile water is contraindicated in patient who is allergic or sensitive to the materials used in the product.
- The use of the product in the patient with active Urinary Tract Infection may exacerbate the condition.

POTENTIAL COMPLICATIONS

- Over-inflation of catheter balloon leading to rupture
- Under-inflation may result in catheter displacement
- Infection risk from compromised sterility
- Discomfort or injury if used improperly

PRECAUTIONS:

- Check the integrity and functionality of the pre-filled syringe before use. Do not use if the package is open or damaged.

INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ PFS with sterile water/ 46 - Rev. 02 Date:11-05-2026

- Handle the device with care & in aseptic condition. Avoid crushing and crimping damage to the product due to application of surgical instruments such as forceps, etc.
- Inspect the product prior to use. Do not use, if the solution is not colourless or contains particles.
- Always use the volume recommended in the Foley catheter instructions for use for inflating the balloon.
- Always ensure the urine is drained from the bladder before the Foley catheter insertion and balloon inflation.
- Prior to balloon inflation using the pre-filled syringe, always ensure that the Foley catheter is fully inserted into the bladder.
- Do not Re-sterilize. Do not Re-use.
- Keep out of reach and sight of children.
- Do use after the Expiry Date.

CAUTION: ⚠

- Carefully read all instructions prior to use.
- Do not use if the package is open or damaged.
- For single use only. Do not Re-use. Discard any partially used product.
- Do not use the product if the pre-filled solution turns hazy, cloudy, or contains any particulated matter.
- The pre-filled solution is sterile and non-pyrogenic if the tip cap is in place.
- Use proper aseptic technique to avoid any risk of contamination.
- Do not use if the pre-filled solution is leaking.

EXPIRATION PERIOD:

Pre-filled Syringe with sterile water has a shelf life of 3 years from the date of manufacturing. The product should not be used after the indicated sterility expiration date.

STORAGE:

Store at controlled temperature at 10 to 40°C (50 to 104°F) and do not freeze.

INSTRUCTIONS FOR USE:

To ensure a sterile application of the product, the health-care worker must comply with the procedures provided for by the health-care institution, applying an aseptic technique to the entire procedure.

- **Step 1.** Carefully inspect the unit packaging. Do not use if the unit packaging is damaged, not intact, or the product is expired.
- **Step 2.** Aseptically peel-open (paper pouch) or tear along the perforation (ribbon pouch) and take out the syringe. Visually inspect the prefilled syringe and its content for clarity,



INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ PFS with sterile water/ 46 - Rev. 02 Date:11-05-2026

discoloration or leakage prior to use. Do not use if any of the above is observed within the syringe.

- **Step 3.** Carefully remove the protective cap and hold the tip of syringe upward and gently depress the plunger to expel air (if any).
- **Step 4.** Aseptically attach the luer connector of the syringe tip to the balloon inflation port of the Foley catheter. (Avoid any contamination)
- **Step 5.** Slowly inject the suggested volume of the solution for accurate and effective anchoring of the balloon in the Foley catheter (**Note:** Always use the volume advised in the Foley catheter instructions for use for inflating the balloon)
- **Step 6.** After use, dispose the empty syringe or any unused in accordance with local or regional laws for medical waste.

INTENDED PATIENT POPULATION:

The Pre-filled Syringe with sterile water is intended to be used in patients irrespective of age and gender.

INTENDED USER:

The Pre-filled Syringe with sterile water is intended to be used by the trained healthcare professionals like Doctors, paramedics, or nursing staff experienced in urinary catheterization.

HOW SUPPLIED:

The Pre-filled Syringe with sterile water is supplied in individually packed in paper pouch / ribbon pouch complete with clear instructions for use and accessories (if any).

WARNINGS: (If space available)

- **For Professional Use Only:** This device must be used by trained healthcare professionals only
- **Not for Injection:** The product is not intended for injection.
- Do not use for drugs reconstitution/dilution, Injection or Rehydration.
- **Sterile Packaging:** Do not use if the packaging is damaged or opened prior to use.
- **Single Use Only:** The pre-filled syringe is designed for single use. Do not reuse or re-sterilize.
- **Proper Handling:** Handle with care to avoid damage to the device, which could compromise its sterility and functionality.
- **Inspect Before Use:** Ensure all components are present and undamaged before use.
- **Follow Instructions:** Carefully follow the included instructions for use to avoid complications
- **Dispose Properly:** Dispose of the device in accordance with local regulations for medical waste.

INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ PFS with sterile water/ 46 - Rev. 02 Date:11-05-2026

RETURN OF DAMAGED PRODUCT

Return the product in its original packing identified by the batch number, purchase information, your reference, and reason for return. Please contact your local distributor office regarding product return/exchange.

DISPOSAL OF THE USED DEVICE

- Disposal of the used product should be done according to hospital procedures, State Regulations, national legislation or regional legislation.
- Device should be disposed after use on individual patient during a single procedure.
- Disposal of the packaging material should be done as biomedical waste.
- Any open packaging should be disposed as biomedical waste.
- Leftover products should be disposed as biomedical waste

REPORTING

In case of any adverse events and/or potentially sight-threatening complications that may reasonably be regarded as product-related and that were not previously expected in nature, severity or occurrence, contact the manufacturer or authorized representative and competent authority of the member state where the user and/or patient is established.

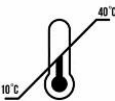
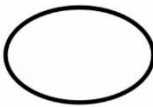
<u>LEGAL MANUFACTURER AND MANUFACTURING FACILITY</u>	<u>AUTHORIZED REPRESENTATIVE</u>		
 Manufacturer: SPM Medicare Pvt, Ltd B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India. Manufacturing Site 1- B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India. Manufacturing Site 2- Industrial Building Plot no 26, Sector 158, Noida, Gautam buddha, Nagar, Uttar Pradesh, 201305. Customer care: No.: +91-840-7070-718 Email: info@spmmedicare.com Web Site: www.spmmedicare.com	<table border="1" style="display: inline-table; vertical-align: top;"> <tr> <td style="padding: 2px;">EC</td> <td style="padding: 2px;">REP</td> </tr> </table> CMC MEDICAL DEVICES & DRUGS S.L. C/ Horacio Lengo, 18 29006, Málaga, Spain Customer Care No.: +34 951 214 054 Email: info@cmcmmedicaldevices.com	EC	REP
EC	REP		

EXPLANATION OF SYMBOLS:

Symbol	Explanation	Symbol	Explanation
	Caution		Catalogue Number

INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ PFS with sterile water/ 46 - Rev. 02 Date:11-05-2026

	Manufacturer		Batch code
	Date & Country of manufacturing		Use-by date
	Do not re-sterilize		Temperature Limit
	Consult instructions for use		Do not use if package is damaged and consult <i>instructions for use</i>
	Keep away from sunlight		Non-Pyrogenic
	Do not re-use		Sterilized using Gamma Irradiation
	Keep dry		Handle with care
	This side up		Maximum stacking limit
	Do not use hook		Recycling
	Medical Device		Model number
	Unique Device Identifier		Single Sterile Barrier System
	Single Sterile Barrier System with protective packaging outside		Biological risks

Revision History:



INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ PFS with sterile water/ 46 - Rev. 02 **Date:**11-05-2026

Rev. No	Reason for Change	Date	Approval
02	Change of Authorized Representative.	11/05/2026	Managing Director