

## **INSTRUCTIONS FOR USE**

**Doc. No.:** SPM/TF- CE/ IVC- Closed/ IFU/ 02- Rev. 04. **Date:** 11/05/2026



---

### **INSTRUCTION FOR USE**

---

#### **CLOSED IV CANNULA**

|                                                                                                         |                                                                                                         |                                                                                                           |
|---------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------|
| <br>DATE: 11-05-2026 | <br>DATE: 11-05-2026 | <br>DATE: 11-05-2026 |
| <b>Prepared By</b><br>Executive- QARA                                                                   | <b>Reviewed By</b><br>Head- RA                                                                          | <b>Approved By</b><br>Managing Director/ Auth. Signatory                                                  |

## **INSTRUCTIONS FOR USE**

**Doc. No.:** SPM/TF- CE/ IVC- Closed/ IFU/ 02- **Rev. 04.** **Date:** 11/05/2026

Follow manufacturer' s instructions for use and institutional procedures for IV Cannula.

**PRODUCT NAME:** IV Cannula (CLOSED IV CANNULA WITH WINGS, WITH EXTENSION LINE & SAFETY FEATURE)

**DESCRIPTION:**

The SPM Closed IV Cannula is a sterile, single-use, peripheral intravenous catheter intended to provide short-term vascular access for the administration of fluids, medications, or blood products.

The device consists of a stainless-steel introducer needle and a flexible PUR (polyurethane) radiopaque catheter designed to provide optimal flexibility, biocompatibility, kink resistance, and clear visibility under imaging. The smooth, bevelled needle enables controlled venous insertion.

The system incorporates a closed infusion mode with an integrated valve that helps prevent blood backflow, leakage, and contamination, thereby reducing blood exposure. A transparent flashback chamber allows quick and accurate confirmation of venous puncture during insertion.

An active safety system is integrated to reduce the risk of accidental needlestick injuries following needle withdrawal, enhancing user safety. Upon confirmation of venous access, the catheter is advanced into the vein and the needle is withdrawn, automatically activating the safety mechanism for safe disposal.

The device includes a stabilizing platform to provide secure fixation and improved handling during insertion. An integrated sliding clamp on the extension line enables controlled flow regulation and helps prevent air entry and backflow.

A Y-connector with two needle-free connectors supports safe multi-line access while reducing the risk of contamination. A removable needle-free valve allows convenient access while maintaining the integrity of the closed system.

The device is color-coded according to gauge size for easy identification. It is manufactured from non-toxic, medical-grade materials and is latex-free.

The Closed IV Cannula is sterilized using Ethylene Oxide (EO), supplied in a validated sterile barrier system, and intended for single use only to support effective infection control and patient safety.

**IMAGE OF THE DEVICE:**



## **INSTRUCTIONS FOR USE**

**Doc. No.:** SPM/TF- CE/ IVC- Closed/ IFU/ 02- **Rev. 04.** **Date:** 11/05/2026

### **DEVICE VARIANTS:**

| <b>Sr. No.</b> | <b>Model</b>                                                                                                             | <b>Product Code</b> |
|----------------|--------------------------------------------------------------------------------------------------------------------------|---------------------|
| 1.             | EASY VEIN CONTROL WAY CLOSED IV CANNULA WITH WINGS WITH EXTENSION LINE & SAFETY FEATURES - PUR CATHETER 14G, With Needle | 3AQ314              |
| 2.             | EASY VEIN CONTROL WAY CLOSED IV CANNULA WITH WINGS WITH EXTENSION LINE & SAFETY FEATURES- PUR CATHETER 16G, With Needle  | 3AQ316              |
| 3.             | EASY VEIN CONTROL WAY CLOSED IV CANNULA WITH WINGS WITH EXTENSION LINE & SAFETY FEATURES- PUR CATHETER 18G, With Needle  | 3AQ318              |
| 4.             | EASY VEIN CONTROL WAY CLOSED IV CANNULA WITH WINGS WITH EXTENSION LINE & SAFETY FEATURES- PUR CATHETER 20G, With Needle  | 3AQ320              |
| 5.             | EASY VEIN CONTROL WAY CLOSED IV CANNULA WITH WINGS WITH EXTENSION LINE & SAFETY FEATURES- PUR CATHETER 22G, With Needle  | 3AQ322              |
| 6.             | EASY VEIN CONTROL WAY CLOSED IV CANNULA WITH WINGS WITH EXTENSION LINE & SAFETY FEATURES- PUR CATHETER 24G, With Needle  | 3AQ324              |
| 7.             | EASY VEIN CONTROL WAY CLOSED IV CANNULA WITH WINGS WITH EXTENSION LINE & SAFETY FEATURES- PUR CATHETER 24G, With Needle  | 3AQ326              |

- ✓ Available in different gauge size with identical colour coded body.

### **INTENDED USE:**

The IV Cannula is a sterile, single use medical device intended to provide peripheral vascular access for short term use.

### **INDICATIONS FOR USE:**

- Phlebotomy (blood sampling)
- Infusion of IV solution
- Provide nutritional support
- Maintain hydration and/or correct dehydration
- Administration of intravenous medications.

### **INTENDED PATIENT POPULATION:**

IV Cannula is intended to be used in patients, irrespective of age and gender. Selection of appropriate gauge size and insertion site shall be determined by the trained healthcare professional based on the patient's age, vein condition, clinical status, and prescribed therapy.

### **INTENDED USER:**

IV Cannula is intended to be used by the trained healthcare professionals like doctors, paramedics, or nursing staff.

## **INSTRUCTIONS FOR USE**

Doc. No.: SPM/TF- CE/ IVC- Closed/ IFU/ 02- Rev. 04. Date: 11/05/2026

### **INTENDED USE ENVIRONMENT:**

IV Cannula is intended to be used in healthcare environments, including hospitals, clinics, and emergency care settings.

### **CLINICAL BENEFITS:**

- Provides reliable short-term peripheral venous access.
- Enables controlled administration of intravenous fluids, medications, blood products, and nutritional solutions as prescribed.
- Reduces repeated venipuncture when ongoing therapy is required.
- Supports maintenance of hydration and therapeutic management through secure venous access.
- The clinical benefits are consistent with the intended purpose and established clinical practice for peripheral intravenous catheterization.

### **PRINCIPLE OF OPERATION:**

The IV Cannula provides access to the peripheral venous system. The introducer needle enables insertion into the vein. After venous entry is confirmed by blood flashback, the flexible catheter is advanced into the vein and the needle is withdrawn.

The catheter hub incorporates a Luer connection compatible with standard infusion systems in accordance with EN ISO 80369-7.

### **PERFORMANCE CHARACTERISTICS:**

The IV Cannula is designed and tested in accordance with EN ISO 10555-1 (Intravascular Catheters – Sterile and Single-Use Catheters – General Requirements) and EN ISO 80369-7 (Luer Connectors).

Performance characteristics include:

- Catheter tensile strength and bond strength
- Leak resistance
- Flow rate performance
- Radiopacity (where applicable)
- Needle integrity and sharpness
- Flashback visualization
- Luer connection compatibility

The device complies with applicable mechanical, physical, and functional performance requirements under the above standards.

## **INSTRUCTIONS FOR USE**

**Doc. No.: SPM/TF- CE/ IVC- Closed/ IFU/ 02- Rev. 04. Date: 11/05/2026**

### **MATERIAL USED:**

| Element                       | Material                             | Composition / Grade                  |
|-------------------------------|--------------------------------------|--------------------------------------|
| Threaded Stopper (Luer Lock)  | High-Density Polyethylene (HDPE)     | Medical Grade HDPE                   |
| Flashback Chamber (Hub Cover) | Polypropylene (PP)                   | Medical Grade PP (110)               |
| Needle Hub                    | Polypropylene (PP)                   | Medical Grade PP (110)               |
| Port Cap (If required)        | Polypropylene (PP)                   | Medical Grade PP (2120)              |
| Silicone Tube                 | Silicone Rubber                      | Medical Grade Silicone Elastomer     |
| Wing Housing (Wing Body)      | Polypropylene (PP)                   | Medical Grade PP (2120)              |
| Needle Cover                  | Polypropylene (PP)                   | Medical Grade PP                     |
| Teflon Holder (Slip Ring)     | Poly Acetal (POM)                    | Medical Grade POM (Polyoxymethylene) |
| Catheter                      | Fluorinated Ethylene Propylene (FEP) | Medical Grade FEP                    |
| Catheter                      | Polyurethane (PUR)                   | Medical Grade PUR                    |
| Needle (Cannula)              | Stainless Steel (SS)                 | AISI 304                             |
| Metal Clip                    | Stainless Steel (SS)                 | AISI 304                             |

- This device does not contain natural rubber latex.
- This device does not contain phthalates, including DEHP.
- The stainless-steel needle may contain trace amounts of nickel. Patients with known hypersensitivity to nickel should be evaluated prior to use.

### **CONTRAINDICATIONS:**

- Hypersensitivity to any material used in the device.
- Not intended for central venous access.
- Avoid insertion in areas of infection, inflammation, oedema, thrombosis, or previously damaged veins.
- Not suitable for administration of highly viscous fluids incompatible with peripheral venous access.

### **POTENTIAL COMPLICATIONS:**

Possible complications associated with peripheral intravenous catheterization may include:

## **INSTRUCTIONS FOR USE**

Doc. No.: SPM/TF- CE/ IVC- Closed/ IFU/ 02- Rev. 04. Date: 11/05/2026

- Local infection or systemic infection
- Phlebitis
- Infiltration or extravasation
- Thrombosis
- Catheter embolism due to catheter damage
- Air embolism
- Hematoma
- Needlestick injury

The likelihood of these complications can be reduced through adherence to proper aseptic technique and institutional protocols.

### **RESIDUAL RISKS:**

Despite design controls, material selection, verification testing, and compliance with applicable standards, certain residual risks such as phlebitis, infiltration, extravasation, thrombosis, infection, catheter embolism, air embolism, or needlestick injury may remain during clinical use.

These risks are inherent to peripheral intravenous catheterization procedures and can be minimized through proper aseptic technique, correct insertion method, secure fixation, monitoring of the insertion site, and adherence to institutional protocols.

All identified residual risks have been evaluated through the manufacturer's risk management process in accordance with EN ISO 14971. The overall residual risk is considered acceptable when weighed against the clinical benefits of the device.

### **INSTRUCTIONS FOR USE:**

1. Inspect the product packaging. Do not use it if the packaging and device is damaged or not intact.
2. Ensure Vent Plug is secure & make sure clamp is not engaged.
3. Twist and remove needle cover.
4. Stabilize the vessel and perform cannulation.
5. Pull back Needle Hub in a control & continuous motion.
6. Look for Initial blood return in the catheter tubing.
7. With the pad of your index finger behind the punch tab, push the catheter forward in the vessel.
8. Look for the continuous blood flow up the extension tube.

## **INSTRUCTIONS FOR USE**

Doc. No.: SPM/TF- CE/ IVC- Closed/ IFU/ 02- Rev. 04. Date: 11/05/2026

9. Stabilize the system and pull back until the needle shield release.
10. Engage the clamp and remove the vent plug.
11. Replace with a needle-free connector or directly the tubing of the infusion set.

### **HOW SUPPLIED:**

The IV Cannula is supplied sterile, single-use, in a validated single sterile barrier system (blister pack) with protective outer packaging.

The sterile barrier system has been validated to maintain sterility for the declared shelf life when stored under recommended conditions.

Each package is labelled with the device description, gauge size, batch number, date of manufacture, expiry date, and UDI carrier.

### **WARNINGS:**

- The use of this product is restricted to a trained doctor or a paramedic.
- In the case of paediatric patients and patients with low blood pressure, flush the extension line with sterile saline prior to use to confirm proper flushing.
- The product should be used according to the instructions for use. SPM Medicare disclaim any responsibility for possible consequences resulting from improper use.
- Store in dry & cool place. Do not expose to heat or direct sunlight.
- The product should not be reprocessed.
- Visually inspect and carefully check the product and packaging before use because improper transport and handling may cause structural and / or functional damage to device or packaging.
- The product is non-toxic, sterile & non-pyrogenic till the package has not been opened or damaged.
- Do not clean or resterilise the product after opening the packet.
- The product should be used immediately after opening the packaging.
- Perform routine monitoring & venipuncture site maintenance according to medical norms.
- Do not reinsert partially or completely withdrawn needle.
- Dispose off the device as per applicable local, state & country laws and regulations.
- Discard after single use.

### **PRECAUTIONS:**

- Inspect the sterile package and device integrity before use. Do not use if the package is open or damaged.
- Select appropriate gauge size based on patient condition and prescribed therapy.
- Maintain strict aseptic technique during insertion and handling.
- Avoid excessive force, bending, crushing, or crimping of the catheter during use.
- Monitor the insertion site regularly for signs of complications.

## **INSTRUCTIONS FOR USE**

Doc. No.: SPM/TF- CE/ IVC- Closed/ IFU/ 02- Rev. 04. Date: 11/05/2026

- Single use only. Do not re-sterilize.
- Do not use after the expiry date indicated on the packaging.

### **CAUTIONS:**

- Store in a cool and dry place
- Discard after single use
- Do not use if the package is open or damaged
- Do not store at extreme temperature and humidity
- Dispose the blood/medicine contaminated product as per local applicable laws.

### **STORAGE:**

- Store at a temperature between 10°C and 40°C in a dry environment.
- Keep away from direct sunlight, moisture, and excessive humidity.
- Do not expose to extreme heat or freezing conditions.

### **STERILITY INFORMATION:**

- This device is supplied sterile and is sterilized using Ethylene Oxide (EO).
- The Sterility Assurance Level (SAL) is  $10^{-6}$ .
- The sterile barrier system has been validated in accordance with EN ISO 11607-1 and EN ISO 11607-2 to ensure maintenance of sterility until the stated expiry date when stored under the recommended storage conditions.
- Do not use if the sterile barrier is damaged or compromised.

### **EXPIRATION PERIOD:**

The IV Cannula has a validated shelf life of **5 years from the date of manufacture**, when stored under the recommended storage conditions stated above.

Do not use the device after the expiry date indicated on the packaging.

The expiry date refers to the sterility and performance of the unopened sterile barrier system.

### **RETURN OF DAMAGED PRODUCT:**

Return the product in its original packing identified by the batch number, purchase information, your reference and reason for return. Please contact your local distributor office regarding product return/exchange.

### **DISPOSAL OF THE USED DEVICE:**

- Disposal of the used product should be done according to hospital procedures, State Regulations, National legislation or regional legislation.
- Device should be disposed after use on individual patient during a single procedure.
- Disposal of the packaging material should be done as biomedical waste.
- Any open packaging should be done as biomedical waste.



## INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ IVC- Closed/ IFU/ 02- Rev. 04. Date: 11/05/2026

- Leftover products should be disposed as biomedical waste.

### REPORTING OF SERIOUS INCIDENTS:

Any serious incident that has occurred in relation to this device must be reported without delay to the manufacturer and to the competent authority of the Member State in which the user and/or patient is established, in accordance with Article 87 of Regulation (EU) 2017/745.

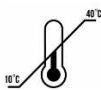
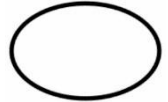
| LEGAL MANUFACTURER AND MANUFACTURING FACILITY                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                          | AUTHORIZED REPRESENTATIVE                                                                                                                                                                                                                                                                                                                                                                                                                                        |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  Manufacturer:<br>SPM Medicare Pvt, Ltd<br>B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India.<br>Manufacturing Site 1- B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India.<br>Manufacturing Site 2- Industrial Building Plot no 26, Sector 158, Noida, Gautam buddha, Nagar, Uttar Pradesh, 201305.<br>Customer care: No.: +91-840-7070-718<br>Email: <a href="mailto:info@spmmedicare.com">info@spmmedicare.com</a><br>Web Site: <a href="http://www.spmmedicare.com">www.spmmedicare.com</a> | <div style="display: flex; align-items: center;"> <div style="border: 1px solid black; padding: 2px; margin-right: 5px;">EC</div> <div style="border: 1px solid black; padding: 2px; margin-right: 5px;">REP</div> <div> CMC MEDICAL DEVICES &amp; DRUGS S.L.<br/> C/ Horacio Lengo, 18<br/> 29006, Málaga, Spain<br/> Customer Care No.: +34 951 214 054<br/> Email: <a href="mailto:info@cmcmmedicaldevices.com">info@cmcmmedicaldevices.com</a> </div> </div> |

### EXPLANATION OF SYMBOLS:

| Symbol                                                                              | Explanation                    | Symbol                                                                              | Explanation      |
|-------------------------------------------------------------------------------------|--------------------------------|-------------------------------------------------------------------------------------|------------------|
|  | Caution                        |  | Catalogue Number |
|  | Manufacturer                   |  | Batch code       |
|  | Date & Country of manufacturer |  | Use-by date      |

## INSTRUCTIONS FOR USE

Doc. No.: SPM/TF- CE/ IVC- Closed/ IFU/ 02- Rev. 04. Date: 11/05/2026

|                                                                                     |                              |                                                                                      |                                                                          |
|-------------------------------------------------------------------------------------|------------------------------|--------------------------------------------------------------------------------------|--------------------------------------------------------------------------|
|    | Do not re-sterilize          |     | Temperature Limit                                                        |
|    | Consult instructions for use |     | Do not use if package is damaged and consult <i>instructions for use</i> |
|    | Keep away from sunlight      |     | Non-Pyrogenic                                                            |
|    | Do not re-use                |    | Sterilized using Ethylene Oxide                                          |
|    | Keep dry                     |     | Handle with care                                                         |
|   | This side up                 |    | Maximum stacking                                                         |
|  | Do not use hook              |  | Recycling                                                                |
|  | Medical Device               |  | Model number                                                             |
|  | Unique Device Identifier     |  | Single Sterile Barrier System                                            |

### Revision History:

| Rev. No | Reason for Change                                                                       | Date       | Approval          |
|---------|-----------------------------------------------------------------------------------------|------------|-------------------|
| 00      | Initial Release                                                                         | 13-04-2024 | Managing Director |
| 01      | Update in instruction                                                                   | 10/26/2024 | Managing Director |
| 02      | MDD to MDR Update                                                                       | 26/05/2025 | Managing Director |
| 03      | Updated Instruction, intended use, Performance Characteristic and Sterility information | 27/02/2026 | Managing Director |
| 04      | Change of Authorized Representative                                                     | 11/05/2026 | Managing Director |